

Visalis™ 100


Reliability in
cataract surgery



Visalis 100

Great value in cataract surgery

With Visalis 100, Carl Zeiss Meditec introduces a durable, predictable and safe phaco system for cataract surgery. Easy, straightforward operation combined with superior ergonomics make Visalis 100 a highly convenient tool to work with. As an integral element of the ZEISS value chain for cataract surgery, Visalis 100 offers clinicians solid returns in reliability and efficiency for predictable surgical and visual outcomes.

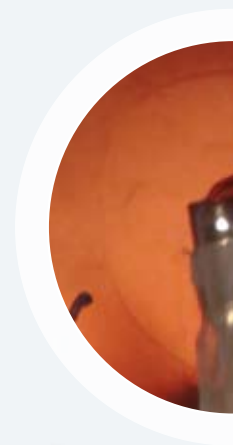



Supporting you every step of the way

Visalis 100 offers excellent usability. Its intuitive graphic user interface with large display, touch-driven control and voice confirmation lets you easily adjust the system. Built for superior handling, Visalis 100 fully supports MICS techniques.

Taking care of business

Thanks to its reliability and stability, the system's servicing and maintenance requirements are very low. What's more, with its compact, lightweight design and small footprint, Visalis 100 requires very little space – qualities that are quite convenient for restricted space in the OR and when needing to transport the system to multiple locations.



Safety as a top priority

Designed for maximum patient safety, Visalis 100 features a completely closed design. There is no internal tubing requiring cleaning, thus eliminating the risk of cross-contamination while increasing the predictability of outcomes. The system's Adaptive Power Control (APC) provides additional safety by delivering surgeons the minimal power necessary during the procedure. The stroke is measured and stabilized by the microprocessor.



*Foot control panel functions
are clearly structured*

ZEISS quality in cataract surgery

*OPMI Lumera® i red reflex technology meets
the phaco efficiency of Visalis 100*

Lens removal with methodical accuracy

Visalis 100 is a highly reliable and stable phaco system, allowing cataract specialists to operate with thorough efficiency and meticulous precision. Its solid performance is the result of the system's premium quality design, including:

- Adaptive Power Control (APC) for predictable surgical and visual outcomes
- Surge Security System to stabilize the anterior chamber
- Fully closed system with no internal tubings for maximum sterility



Technical data

Visalis 100

Dimensions

Height: 193 mm, width: 400 mm, depth: 405 mm

Weight

13.0 kg

Irrigation

- Gravity fed fluid delivery; eye pressure determined by height of irrigation source
- Solenoid valving element

Suction

- Peristaltic aspiration pump
- 5 to 500 mg HG vacuum range
- Flow rate from 2 to 50 cc/min
- 0.5 to 12 s aspiration rise time

Phacoemulsification

- Piezoelectric handpiece with approx. 40 kHz
- Panel or linear control of U/S power with system footswitch

Vitrectomy

- Adjustable cutting speed from 60 to 700 cuts per minute
- Single cut mode available

Diathermy

- Adjustable bipolar power from 5 to 100%
- Linear control with system footswitch



CE 0297

Carl Zeiss Surgical GmbH

A Carl Zeiss Meditec Company
73446 Oberkochen
Germany

Fax: +49 (0) 73 64/20-48 23
surgical@meditec.zeiss.com
www.meditec.zeiss.com/eye

EN_30_010_11011 Printed in Germany AW-CZ-V/2009 Too
Visalis is a trademark of Carl Zeiss. The contents of the brochure may differ from the current status of approval of the product in your country. Not for sale in the US.
Please contact our regional representative for more information. Subject to change in design and scope of delivery and as a result of ongoing technical development.
Printed on elemental chlorine-free bleached paper.
© 2009 by Carl Zeiss Meditec AG. All copyrights reserved.